ITSU

User Manual
IS0131 使用説明書

ITSU

User Manual

All information is correct at the time of print

CONTENTS

Safety Instructions	03
Warranty Claims	05
Main Unit & Accessories	06
Exchangeable heads	06
Instructions	07
Product Maintenance	08
Environmental safety	08
Troubleshooting	09
Product Specifications	09
Contact us	10

Safety Instructions

i. Usage environment

- Do not use the appliance in a wet condition or dusty environment as this may result in malfunction or electric shock.
- It can only be used on dry skin and cleansed skin surface. It should not be used across clothing. Gently press and move on the skin for approximately 60 seconds per part.
- Use the device only on the soft tissues of the body without causing pain or discomfort. Do not use it on the head or any hard or boned part of the body.
- Use only the vibrating head that produces the best results.
- Whether it is controlled pressure or pressure, it can produce bruising. Always check the physiotherapy area but signs of pain or discomfort should be discontinued immediately.
- Keep your fingers, hair or other parts of the body away from the shaft and back of the vibrating head, as there is a possibility of being pinched.
- Place the appliance in a safe location.
- Do not expose the appliance to a heater or direct sunlight.
- This product is only intended for personal use only. Commercial use is prohibited.

ii. Health precautions

If you are experiencing health complications or conditions as stated below, please consult your doctor before use:

- Suffering from heart disease.
- Pregnancy.
- Any health concerns.

iii. When the product should not be used

- If the use exceeds 25 minutes at one time.
- This unit should not be used by children or invalids without adult supervision.
- Before sleeping.
- If water is accidentally spilled on to the appliance.
- If the user detects any malfunction or other abnormal conditions during operation.

iv. Caution and warning

- Never let the device run or charge unattended.
- Do not use any attachments that are not recommended by ITSU.
- Do not operate this appliance if it has a damaged cord or plug. Please return the unit to ITSU. Service Center for examination and repairs.

Safety Instructions

- Do not insert your hands or any foreign matters into the massager's vents and gaps.
- Do not drop or improperly use the device.
- Charge only with 16.8 VDC device charger.
- Do not change or modify the device in any way.
- Do not cover the appliance during operation. As this may result in malfunction or fire.
- Do not crush the appliance.
- Do not overload the electrical outlet.
- If the device malfunctions, disconnect it from the power source and contact customers services.
- Do not let the device come in contact with metal object, such as belt buckle or jewelry.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.

V. Note

- To ensure safe and proper use of the appliance, do not operate it before reading this user manual carefully first.
- It is normal to hear the noises during use, because they are resulted from vibration of product.
- If the appliance is stored in a cold place and being taken into a warm environment please wait one hour before use to get the optimum massaging level.
- The appliance is only to be used with power accessories provided.
- Cleaning and user maintenance shall not be made by children without supervision.

Warranty Claims

1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from ITSU.
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
6. Warranty does not cover accessories and attachments which do not belong to the massage unit.
7. Warranty is not effective to rental, business, commercial, institutional or other non-residential users.
8. All services covered by this warranty must be approved by ITSU and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, ITSU reserves the right to make substitution in lieu of repair or replacement.


Main Unit & Accessories

- a) Handheld Massager
- b) 4 types of exchangeable massage heads
- c) USB Charging cable

Handheld massager

- 1. ON/OFF switch 
- 2. Detachable, replaceable shock head.
When removing vibrating head, pull it out with ease.
- 3. USB charging cable .
- 4. Timer setting is 15 minutes.

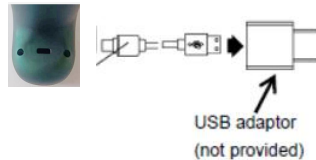Control panel

Power Key:

-Long press for 2 seconds to ON/OFF the massager.
short press to adjust massage level or resume from pause condition.

Power Indicator:

- Signal light:. The signal light will stay illuminated once it is turn ON.
-The four indicators light up once in sequence.
-When charging. The signal lights flash in turn, and the four lights are fully illuminated once fully charged.
- When the power is insufficient, the first power indicator lamp flashes, please charge it immediately.

Level Indicator

Level Indicator:

- (1) Level 1 , suitable for warming up muscles.
- (2) Level 2 , suitable for relaxing the fascicles.
- (3) Level 3 , suitable for lactic acid removal.
- (4) Level 4 , suitable for deep massage.
- (5) Level 5 , suitable for professional massage.
- (6) Level 6, Intelligent speed regulation: Five Gear Lights Rolling Display. Automatic levels conversion according to your pressure, Meaning to say , the strength and speed of massage will be regulated automatically.

Exchangeable heads

U head

Ball head

Flat head

Bullet head

Instructions

Power operation cautions

- Long press the power switch(about 2 second) on the side of the machine to turn it on and start the device (the power light and level light is on).
- Long press again to turn off the device.

Device

- Press  to adjust massage level (speed).Default speed is at level 1.Press again to adjust to level 2,and so on.There are 6 levels.

Charging

- Before the first use, charge it for about 3 hours with 5V 2A power source.
- When the power is insufficient, the first power indicator lamp flashes
- When charging, the corresponding power indicator light flashes.
- After the charge is completed, the four-unit power lamp is always on.
- It is not recommended to completely drain the battery.
- Make sure the power is off when not in use. The operating time of the device depends on the speed level during use and the applied pressure.

Timing function

The device runs for about 15 minutes.

Product Maintenance

TO STORE

- Place the massager in the suitcase provided or in a safe, dry and cool place. Avoid contact with sharp edges or pointed objects that might cut or puncture the surface.

TO CLEAN

- Avoid immersing the unit in any liquid during cleaning.
- Unplug the unit and allow it to cool before cleaning.
- Clean only with a soft, slightly damp sponge.
- Do not allow water or any other liquid to go inside the unit.
- Avoid using abrasive cleaners, brushes, gasoline, kerosene, glass/ furniture polish or paint thinner to clean the product.

CAUTION

- Avoid tearing on the appliance surface.
- Only gentle force should be exerted against the unit in order to eliminate any risk of injury.

Clarication on environmental safety

This product contains a lithium battery and must not be disposed of with normal household waste and must be disposed of safely at a designated hazardous waste or electronic product collection point.

Troubleshooting

This massager is controlled by an integrated chip and is susceptible to power surges from +10% to -15%.Due to this, the massager might display the following:

Conditions	Solutions
Not operating at all	1.Please rechargeable before first time use. 2.Please swich ON the product.
Speed too high or too low	1.Press Mode to adjust massage intensity. 2.There are 6levels.
Product malfunction	1. Let product to cool down for 30 minutes and turn ON again . If still fail,please contact nearest service centre..

Product Specifications

Model of product	IS0131
Power supply	Type-c usb charging
Power consumption	5V 2A
Dimensions (length/ width/ height)	147x126x50mm
Weight	380g

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ITSU World(HK) Limited