ITSU

User Manual
IS2005 使用説明書

ITSU

User Manual
All information is correct at the time of print

CONTENTS

Safety Instructions	03
Warranty Claims	05
Main Unit & Accessories	06
Installation Guid	07
Installation steps	08
Product Operation Guide	09
Product Maintenance	11
Troubleshooting	12
Product Specifications	13

Safety Instructions

i. Usage environment

- Do not use the appliance in a wet condition or dusty environment as this may result in malfunction or electric shock.
- Place the appliance in a safe location.
- Do not expose the appliance to a heater or direct sunlight.
- This product is only intended for personal use only. Commercial use is prohibited.

ii. Health precautions

If you are experiencing health complications or conditions as stated below, please consult your doctor before use:

- Suffering from heart disease.
- Pregnancy.
- Any health concerns.

iii. When the product should not be used

- If the use exceeds 15 minutes at one time.
- This unit should not be used by children or invalids without adult supervision.
- Before sleeping.
- If the fabric cover is torn.
- If water is accidentally spilled onto the appliance.
- If the user feels abnormal or severe pain during the massage.
- If the user detects any malfunction or other abnormal conditions during operation.
- If there is lightning.

iv. Caution and warning

- Unplug the appliance from the outlet when not in use.
- Do not use any attachments that are not recommended by ITSU.
- Do not operate this appliance if it has a damaged cord or plug. Please return the unit to ITSU. Service Center for examination and repairs.
- Do not insert your hands or any foreign matters into the massager's vents and gaps.
- Handle heated surfaces carefully. The unattended use of the heated surfaces by children or invalids may be dangerous.
- This product is only intended for personal use only. Commercial use is prohibited.

Safety Instructions

- Do not cover the appliance during operation. As this may result in malfunction or fire.
- Do not crush the appliance and avoid any sharp folds.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

V. Note

- This appliance can be used by children at the age of 8 or above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under guidance or if they have been informed of the risks involved and instructions concerning use of the appliance in a safe way.
- Children shall not play on the appliance.
- Cleaning and maintenance shall not be performed by children without supervision.
- To ensure safe and proper use of the appliance, do not operate it before reading this user manual carefully first.
- It is normal to hear the noises during use, because they are resulted from vibration of product.
- For people under medication or with a medical condition, please consult your doctor before use. Please do not use the product if you are not feeling well.
- If the appliance is used in a cold room, do not increase the room temperature abruptly. It is advised that the temperature is gradually increased to a normal level.
- If the appliance is stored in a cold place and being taken into a warm environment, please wait one hour before use to get the optimum massaging level.
- Do not remove the cover of the appliance. Mishandling the internal components may result in malfunction or electric shock.
- Do not connect or disconnect the power plug from the electrical outlet with wet hands.
- Do not use the appliance if the electrical outlet is loose.
- The appliance is only to be used with power supply provided.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.



Warranty Claims

1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from ITSU .
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
6. Warranty does not cover accessories and attachments which do not belong to the massage unit.
7. Warranty is not effective to rental, business, commercial, institutional or other non-residential users.
8. All services covered by this warranty must be approved by ITSU and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, ITSU reserves the right to make substitution in lieu of repair or replacement.

Main Unit & Accessories

- A) Massager
- B) Power adapter



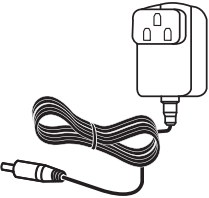
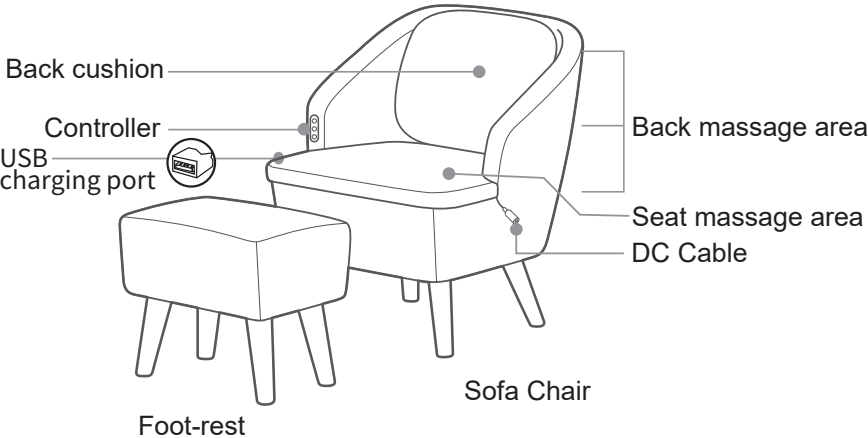
Controller Diagram

Power adapter

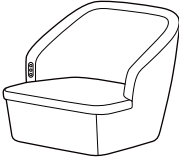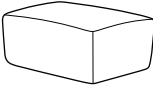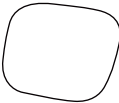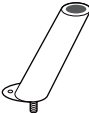
Installation Guid

Parts list:


NO.	Parts name	Pic	Qty
1	Main sofa		1
2	Main foot-rest		1
3	Back cushion		1
4	Solid wood chair leg		8
5	Allen wrench		1
6	Socket Head Cap Screw		9
7	Screw washer		8
8	Chair leg felt Pad		4

Installation steps


1 up side down main sofa

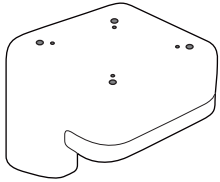
2 put washer & chair leg screw together

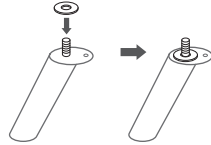
3 rotate clockwise to install chair leg

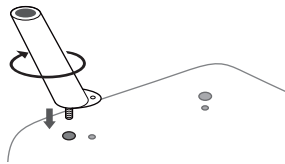
4 lock socket head cap screw into chair with Allen wrench

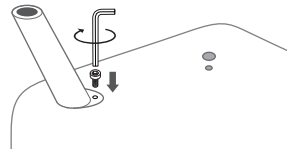
5 install rest of chair legs one by one

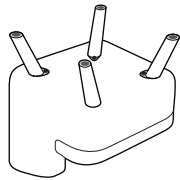
6 sofa chair installation complete

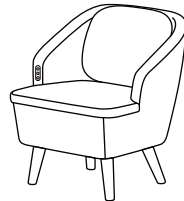7 follow previous steps to install Foot-rest

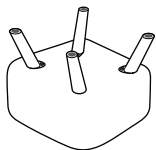8 Foot-rest installation complete

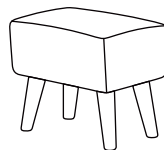Product Operation Guide

i. Product features

- 6 smooth massage balls with kneading and tap to relax waist muscles.
- Physical therapy function and massager stimulation to soothe multiple aches and provide a more thorough massage.
- Versatile control button, handy.
- After using for 15 minutes, the power will automatically end.
- The ergonomic design allows the product to be used in different occasions and places.

ii. Operation instruction

a) Preparation before use

- Check if there is any damage to the appliance.
- Place the appliance on the part which needs to be massaged, before ensuring that the appliance is properly placed.

b) Usage method

1. Connect the DC head of the power adapter to the DC connector of the appliance, and plug the other side into the power base.

- Default timer is 15 minutes .

- Key operation:

Long press 2S for ON/OFF function, default massage mode 1
*press again to switch model 2, third press switch model 1

Press to enable/disable back&seat HEAT function

Press to enable/disable back&seat VIBRATION function

Product Maintenance

TO STORE

- Place the massager in this box or in a safe, dry and cool place. Avoid contact with sharp edges or pointed objects that might cut or puncture the surface.
- Please wind up and store the power cord under the appliance bottom frame. Do not hang the unit by the power cord.

TO CLEAN

- Avoid immersing the unit in any liquid during cleaning.
- Unplug the unit and allow it to cool before cleaning.
- Clean only with a soft, slightly damp sponge.
- Do not allow water or any other liquid to go inside the unit.
- Avoid using abrasive cleaners, brushes, gasoline, kerosene, glass/ furniture polish or paint thinner to clean the product.
- Do not attempt to repair the product. For servicing, please contact ITSU.

CAUTION

- Avoid tearing on the appliance surface.
- Only gentle force should be exerted against the unit in order to eliminate any risk of injury.
- Beware when put on light color PU/fabric seat, some material may have color dyeing
- Beware when put printing document on product, some material may have color dyeing


Troubleshooting

This massager is controlled by an integrated chip and is susceptible to power surges from +10% to -15%. Due to this, the massager might display the following:

Conditions	Possible reason	Solutions
NO Inoperation	Poor power outlet	1.Check the power outlet for electrical output 2.Check whether the outlet plug is in good contact.
Power socket is good but not working	fault	Please unplug the plug and contact the seller for repair.

Product Specifications

Model number of the product	IS2005
Power source	100-240V~
Power frequency	50/60Hz
Power consumption	36W
Voltage	12V= 3A
Timer	15minutes
Dimension (Length/width/height)	603x570x720mm
Weight	13.3kg

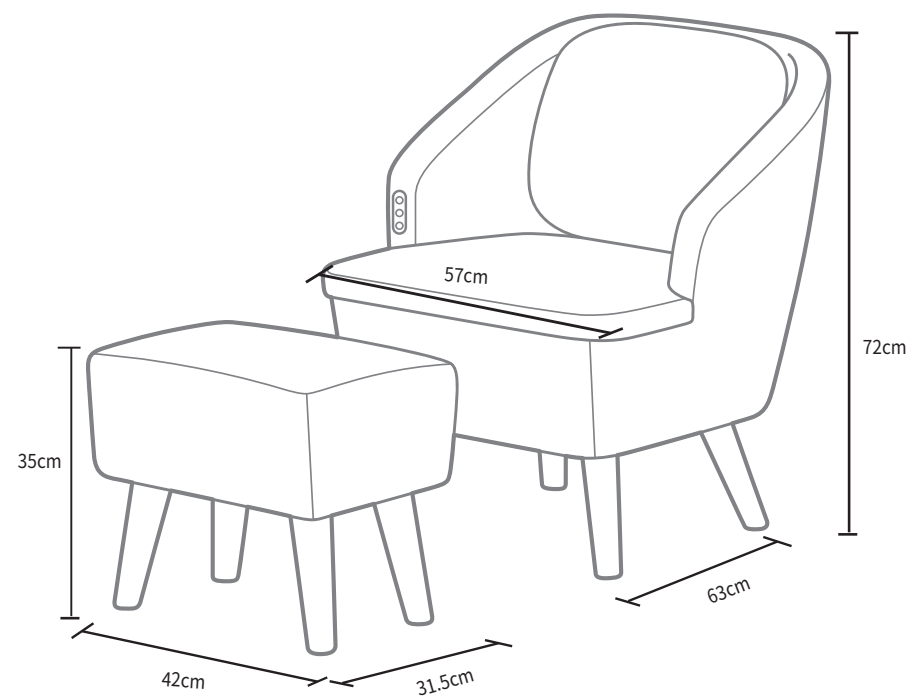Contact Us

ITSU WORLD SDN. BHD.
Unit BO03-B-13A-1,
Menara 3A, KL Eco City,
No.3, Jalan Bangsar,
59200 Kuala Lumpur.
Tel: 1700-81-4878

ITSU World(HK) Limited
Flat A6, 20th Floor, TML Tower, No. 3
Hoi Shing Road, Tsuen Wan, N.T.
Tel: +852 3611 2228

PT Healthy World
Lodan Center Block O No 20
Jalan Lodan Raya No 2
Jakarta Utara 14430
Tel: +62 21-6983 1123

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